

DIVAN ÉLECTRIQUE "LEMI 3" COUCH ELECTRIC "LEMI 3"

Réf : M018210

- Pratique et universel.
- Lit/fauteuil à élévation oléo dynamique, actionné par une pompe hydraulique à pédale de qualité, avec la possibilité de tourner sur 360°. Le basculement (trendelenburg) et les mouvements du dossier, des jambes et de la tête sont réglables manuellement, le dossier est mouvementé par l'intermédiaire d'un vérin à gaz, ce qui permet à l'opérateur de trouver facilement la bonne position même avec des patients âgés, garantissant sécurité et confort.
- La structure en acier électro zingué, totalement fermée par un carter en abs, confère élégance et sécurité.
- La version standard est livrée avec un matelas plat constitué d'une matière synthétique indéformable.

- Practical and universal.
- Bed / armchair with dynamic oleo lift, operated by an hydraulic pedal pump of optimum quality, with possibility of turning 360°. The tilting (trendelenburg) and the movements of the backrest, legs and headrest are manually adjustable, the backrest is moved by of a gas cylinder, allowing the operator to easily find the right position even with the elderly patients, guaranteeing safety and comfort.
- The zinc-plated steel structure, completely closed by an ABS casing, gives elegance and security.
- The standard version comes with a flat mattress made, non-deformable synthetic material

Réf.	Long./Length	Larg./Width	Haut./Heigh
M018210	190 cm	74 cm	86 cm



EN OPTION

OPTIONAL

Pédale

Pedal

Réf. : M018211



Lemi

GAMME COULEURS POUR LE LEMI 3 ET LEMI 4


COULEURS REVÊTEMENTS COATINGS COLORS



6



1155

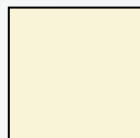

1182



1957



Noir



Gris



Blanc

COULEURS STRUCTURES STRUCTURES COLORS